LOGIKO

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Venduto a _____

Via _____ Località _____

Honsun (Nantong) Co. Ltd, No.8 Tongxing
Rd. Nantong Economic & Technological Dev. area, 226009, Jiangsu Province, China

EC REP Shanghai International Holding Corp. GmbH, Eifestrasse 80, 20537 Hamburg, Germany

Distribuito da / Distributed by:

MORETTI S.P.A.

Via Bruxelles, 3 - Melegnano 20022 Cavriaglia (Arezzo) Tel. +39 055 96 21 11

www.morettispa.com email: info@morettispa.com

MADE in P.R.C.

Class I Medical Device with measuring function
Italian Legislative Decree no. 46 dated 24/02/97 implementing Directive EEC 93/42
and subsequent amendments

CODES

DM344	Palm aneroid sphygmomanometer kit with 3 cuffs
DM349	Standard double tube palm aneroid sphygmomanometer
DM361	Wall aneroid sphygmomanometer with basket
DM362	Wall aneroid sphygmomanometer with tank
DM366	Stand type aneroid sphygmomanometer with basket
DM367	Stand type aneroid sphygmomanometer with tank
DM370	Desk aneroid sphygmomanometer

INTRODUCTION

Thank you for having chosen a **LOGIKO** aneroid sphygmomanometer by Moretti.

The **LOGIKO** sphygmomanometers by Moretti are designed and built to meet all your demands for safe, practical, correct use. This manual provides some suggestions as to how to correctly use the device you have chosen and gives some valuable advice for your safety. Please read through the manual carefully before using the aneroid sphygmomanometer. Contact your retailer directly for questions or further assistance.

INTENDED PURPOSE

The **LOGIKO** aneroid sphygmomanometers by Moretti are devices used to measure blood pressure. They are based on the auscultatory technique according to the Riva Rocci/Korotkoff method.

GENERAL DESCRIPTION

The air pressure within the cuff needed to take the measurement is created using a bulb. The bleed valve allows for fine regulation of the air flow, allowing for measurement.

Each manometer has a graduated scale with a measurement range from 0 to 300 mmHg (mmHg = millimetres of mercury), marked out into 2 mmHg notches and greater notches every 10 mmHg for easy reading. Reading accuracy for the whole range is ± 3 mmHg.

With this system, both the **SYSTOLIC** and **DIASTOLIC** pressure are measured.

WHAT IS BLOOD PRESSURE?

Blood pressure is the pressure exercised by the blood on the blood vessel walls. The maximum arterial pressure during the cardiac cycle is defined as **SYSTOLIC** pressure. The lowest arterial pressure is defined as **DIASTOLIC** pressure.

WHAT IS NORMAL BLOOD PRESSURE?

Blood pressure can vary greatly because it depends on a great many factors. Generally speaking, blood pressure is lower in the summer and higher in the winter. Blood pressure can vary with atmospheric pressure and is significantly affected, for example, by physical demands, emotions, stress, meals, medicines, alcohol, smoke, age, etc..

It is a good idea to write down daily measurements and then consult a doctor in order to suitably define a normal blood pressure for your own body.

N.B. Blood pressure varies with age, you will therefore need to consult a doctor to find out what blood pressure is normal for you. Under no circumstances should you alter the dosage of any drugs your doctor may have prescribed!

REGULATIONS AND DIRECTIVES OF REFERENCE

The **LOGIKO** sphygmomanometers by Moretti are built in compliance with current standard EN 1060-1-2:

EN 1060-1 Non-invasive sphygmomanometers, general requirements;

EN 1060-2 Non-invasive sphygmomanometers, supplementary requirements for mechanical sphygmomanometers.

WARNINGS FOR USE

Please read carefully this manual to ensure correct use of the device. Always follow the instructions given with regard to the product you have purchased. Always keep the device out of the reach of children.

Never take pressure above 300 mmHg.

**GENERAL ASSEMBLY OF PARTS**

DOUBLE TUBE ANEROID SPHYGMOMANOMETER - Connect one tube of the air chamber to the manometer and the other tube to the bulb with the pressure regulation valve.

SINGLE TUBE PALM ANEROID SPHYGMOMANOMETER - Connect the tube of the air chamber to the relevant plug on the palm aneroid sphygmomanometer.

DOUBLE TUBE PALM ANEROID SPHYGMOMANOMETER - Connect one tube of the air chamber to one of the two plugs on the palm sphygmomanometer. Connect the other tube to the remaining plug on the sphygmomanometer.

N.B. To assemble the desk and wall-mounted sphygmomanometer, please refer to the relevant additional instructions provided with this manual.

LOGIKO sphygmomanometers are supplied with standard adult cuffs (arm circumference from 26 to 36 cm). The following cuffs are also available to order:

- Children's cuff (arm circumference from 18 to 26 cm)
- Obese cuff (arm circumference from 36 to 46 cm)



The palm aneroid sphygmomanometer (DM344) is always supplied as a kit containing all 3 cuff types (children, adult and obese).

Before using the sphygmomanometer, you need to choose the right cuff as this will help guarantee greater accuracy in measuring your blood pressure.

Please contact your retailer if in doubt as to which cuff to use.

HOW TO USE

The user should be in a warm place. Clothing must be removed from the left arm. Sit at a table or desk, where it is easier to rest your arm. Position your elbow on the table so that it is approximately the same height as your heart. Turn your arm so that the palm of your hand faces upwards. Wrap the cuff around your bare arm, ensuring that the centre of the cuff corresponds to the brachial artery. The lower end of the cuff must be approximately 2-3 cm above the elbow joint. Tighten the cuff so that it adheres closely to your arm. You must be able to insert a finger between the cuff and your arm.

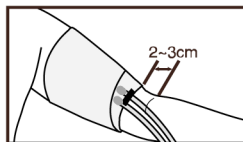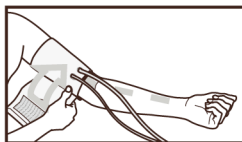
N.B. Do not tighten too far or the blood flow will be reduced significantly.

Insert the phonendoscope head (OPTIONAL) beneath the centre of the cuff, approximately 3-4 cm away from the elbow joint.

N.B. The phonendoscope should not be pressed too firmly against the cuff, or the diastolic pressure measured will be a false reading.

Wear the phonendoscope earpiece in order to hear the heartbeat. Now inflate the air chamber of the cuff using the pump supplied. To do so, first ensure that you have closed the pressure bleed valve. Pump air into the circuit until pressure reaches approximately 30 mmHg above normal systolic pressure (if this is the first time you measure your blood pressure, inflate up to approximately 180 mmHg). Remain very still whilst measuring pressure. Start to deflate the air chamber slowly, adjusting the bleed valve in order to obtain a constant loss of pressure of approximately 2-3 mmHg per second. Listen to the phonendoscope as you watch the manometer. When you can clearly perceive at least two consecutive heartbeats, this is your peak (SYSTOLIC) pressure. When these beats cease being recognisable, you have your lowest (DIASTOLIC) pressure. After having obtained both pressure readings, open the bleed valve completely to quickly empty the air chamber of all residual pressure. Now remove the phonendoscope and cuff.

Do not keep the cuff inflated for too long.

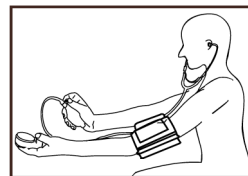WARNINGS FOR MEASURING BLOOD PRESSURE

- Generally speaking, you should always use the left arm.
- Uncover your arm, remove all clothing from it.
- Avoid all activity prior to taking the measurement.
- Rest for 5 to 10 minutes prior to measuring blood pressure.
- Do not talk, eat, drink or move whilst taking your blood pressure.
- Do not press the head of the phonendoscope with your fingers. This could affect the blood pressure reading.
- Always measure blood pressure at the same time.
- Always record the time and date of measurement.

CARE AND MAINTENANCE

The **LOGIKO** devices by Moretti are CE marked and carefully checked before being released for sale.

Each sphygmomanometer is calibrated and subjected to a final thorough control.

Before each use, while resting, check that the needle on the aneroid sphygmomanometer dial is within the tolerance range permitted of ± 3 mmHg (this can be noted by the greater marking of the graduated scale slightly above and below zero). Appliances out of this tolerance range must be recalibrated. Consult your retailer for qualified technical assistance. Only original replacement parts and accessories must be used for repairs.

Pay careful attention to the following:

- Pressure must not exceed 300 mmHg.
- Do not leave the device in damp places.
- Do not drop the device (should it fall, always ensure that an authorised laboratory checks that the device is functioning correctly).
- Do not leave the device in direct sunlight.
- Keep away from all heat sources.

CLEANING

MANOMETER - Only use a soft, dry cloth.

BULB WITH VALVE - Only use a soft damp cloth and then dry thoroughly.

CUFF - To correctly clean the cuff, remove the air chamber, then wash the cuff using warm water and neutral soap.

AIR CHAMBER - Wash with warm water and neutral soap, taking care not to allow the water to enter. Should a little water enter the air chamber, allow to dry well before use by arranging the tubes pointing downwards to allow the water to drain out.

PHONENDOSCOPE (OPTIONAL) - Only use a soft damp cloth and then dry thoroughly.

For the remaining parts, use warm water and neutral soap. Dry well.

GENERAL WARNINGS

Do not use the following product for purposes other than that defined by this manual.

MORETTI S.P.A. refuses all liability for damages caused by improper use of the device or by its use for a purpose other than that specified in this manual.

The manufacturer reserves the right to make changes to the device and this manual in order to improve the device characteristics without advance notice.

TECHNICAL CHARACTERISTICS

Storage temperature:	from -20°C to +70°C
Humidity:	85% max. relative
Operative temperature:	from 0°C to 46°C
Measurement range:	from 0 to 300 mmHg
Minimum scale:	2 mmHg
Measuring accuracy:	± 3 mmHg from 18°C to 33°C ± 6 mmHg from 34°C to 46°C
Air leakage:	< 4 mmHg/min.

SPARE PARTS / ACCESSORIES

Please refer to the general MORETTI catalogue for spare parts and accessories

N.B. The table below provides information on how to correctly use the cuff.

CUFF TYPE

ARM CIRCUMFERENCE IN CM

Children	18 - 26
Adults	26 - 36
Obese	36 - 46

ASSEMBLY INSTRUCTIONS **DM366**

1. Insert the telescopic tube (ref. 7) into the wheeled stand (ref. 8)
2. Connect the two parts using the screw (ref. 10) and seal (ref. 9)
3. Tighten the screw (ref. 10) using the wrench supplied (ref. 11)
4. Insert the manometer (ref. 12) into the upper rod of the telescopic tube, applying pressure in a downwards direction;
5. Adjust manometer inclination using the screw and nut (ref. 4 and ref. 5)
6. Adjust the stand height using the adjustment ring (ref. 6)
7. Connect one end of the coil tube to the plug on the back of the manometer and the other end to one of the cuff tubes
8. Connect the other cuff tube to the bulb with the valve

ASSEMBLY INSTRUCTIONS **DM361**

1. Choose a fixing point (e.g. a wall) that is sufficiently resistant and secure for assembling the manometer
2. Establish the height at which the manometer is to be mounted so as not to constitute a danger and/or hindrance to people or objects
3. Make two holes in the wall to fit the plugs supplied (6 mm)
4. Screw in the manometer bracket using the two screws supplied
5. Mount the manometer onto the bracket using the screw and nut (ref. 4 and ref. 5). Adjust inclination before blocking
6. Connect one end of the coil tube to the plug on the back of the manometer and the other end to one of the cuff tubes
7. Connect the other cuff tube to the bulb with the valve

ASSEMBLY INSTRUCTIONS **DM367**

1. Insert the telescopic tube (ref. 5) into the wheeled stand (ref. 6)
2. Connect the two parts using the screw (ref. 8) and seal (ref. 7)
3. Tighten the screw (ref. 8) using the wrench supplied (ref. 9)
4. Insert the manometer (ref. 10) into the upper rod of the telescopic tube, applying pressure in a downwards direction
5. Adjust manometer inclination using the screw and nut supplied
6. Adjust the stand height using the adjustment ring (ref. 4)
7. Connect one end of the coil tube to the plug on the back of the manometer and the other end to one of the cuff tubes
8. Connect the other cuff tube to the bulb with the valve


ASSEMBLY INSTRUCTIONS **DM362**

1. Choose a fixing point (e.g. a wall) that is sufficiently resistant and secure for mounting the manometer
2. Establish the height at which the manometer is to be mounted so as not to constitute a danger and/or hindrance to people or objects
3. Dismantle the rear bracket of the basket (ref. 2), acting on the internal screw
4. Make two holes in the wall to fit the plugs supplied (6 mm)
5. Screw in the manometer bracket using the two screws supplied
6. Re-assemble the manometer onto the bracket using the internal screw. Adjust angle before blocking
7. Connect one end of the coil tube to the plug on the back of the manometer and the other end to one of the cuff tubes
8. Connect the other cuff tube to the bulb with the valve

DM366 DESK ANEROID
SPHYGMOMANOMETER W/BASKET

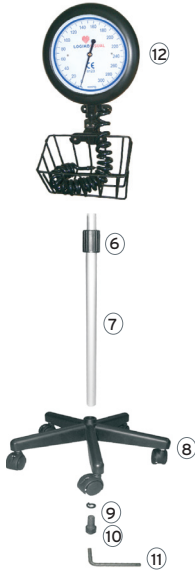DM367 ANEROID SPHYGMOMANOMETER
STAND TYPE W/TANK

DM361 ROUND WALL-MOUNTED ANEROID
SPHYGMOMANOMETER W/BASKET

- 1 Wall-mounted manometer (back)
- 2 Basket
- 3 Coil connection tube
- 4 Screw
- 5 Nut
- 6 Adjustment ring
- 7 Telescopic tube
- 8 5-leg stand with wheels
- 9 Seal
- 10 Screw
- 11 Wrench
- 12 Stand type manometer

- 1 Wall-mounted manometer (back)
- 2 Basket
- 3 Coil connection tube
- 4 Adjustment ring
- 5 Telescopic tube
- 6 5-leg stand with wheels
- 7 Seal
- 8 Screw
- 9 Wrench
- 10 Stand type manometer

ENGLISH

13

WARRANTY

This product is guaranteed for **2 YEARS** as from the date of purchase. The distributor reserves the right to repair or replace any parts affected by manufacturing flaws or materials used, at no expense to the customer. **The warranty shall not apply in the event of damages caused by improper use, abuse, alterations or dismantling by unauthorised persons.** For repair, in compliance with the above conditions for warranty application, please send the product, together with the purchase receipt, to the nearest local retailer.

MORETTI®

LOGIKO

WARRANTY CERTIFICATE

Product _____

Purchased on (date) _____

Retailer _____

Address _____ Town/city _____

Sold to _____

Address _____ Town/city _____

Honsun (Nantong) Co. Ltd, No.8 Tongxing
Rd. Nantong Economic & Technological Dev. area, 226009, Jiangsu Province, China

Shanghai International Holding Corp. GmbH, Eiffestrasse 80, 20537 Hamburg, Germany

Distribuito da / Distributed by:
MORETTI S.P.A.

Via Bruxelles, 3 - Melegnano 52022 Cavriglia (Arezzo) Tel. +39 055 96 21 11
www.morettispa.com email: info@morettispa.com
MADE in P.R.C.

14

Note

[illegible]

Note

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

Distribuito da / Distributed by:
Via Bruxelles, 3 - Melegnano
53022 Casavola (Arezzo)
Tel. +39 055 96 21 11

Fax. +39 055 96 21 200
www.morettipa.com
info@morettipa.com